

## HUMAN SOLUBLE NEURAL CELL ADHESION MOLECULE (NCAM-1/CD56) ELISA KIT

FOR THE QUANTITATIVE DETERMINATION  
OF HUMAN SNCAM-1 CONCENTRATIONS  
IN CSF AND PLASMA



FOR RESEARCH USE ONLY. NOT FOR USE IN  
DIAGNOSTIC PROCEDURES.

### PURCHASE INFORMATION:

ELISA Name	Human Soluble NCAM-1 ELISA
Catalog No.	SK00266-01
Lot No.	
Formulation	96 T
Standard Range	47-3000 pg/mL
Sensitivity	24.5 pg/mL
Sample Volume	100 µl
Sample Type	CSF & EDTA plasma
Dilution Factor	Optimal dilutions should be determined by each laboratory for each application
Specificity	Human NCAM-1 only
Intra-assay Precision	4-6%
Inter-assay Precision	8-10%
Storage	2°C-8°C

### Order Contact:

**AVISCIERA BIOSCIENCE INC.**

**2348 Walsh Ave., Suite C**

**Santa Clara, CA 95051**

**Tel: (408) 982 0300**

**Fax: (408) 982 0301**

**Email: [Sales@AvisceraBioscience.com](mailto:Sales@AvisceraBioscience.com)**

**[Info@AvisceraBioscience.com](mailto:Info@AvisceraBioscience.com)**

**[www.AvisceraBioscience.com](http://www.AvisceraBioscience.com)**

**INTRODUCTION**

Human soluble NCAM-1 immunoassay is a 3.5 - 4.5 hour solid phase ELISA designed to measure human NCAM-1 in CSF and plasma. It contains recombinant human NCAM-1 and antibodies raised against this protein. It has been shown to accurately quantify recombinant human NCAM-1. Results obtained with naturally occurring NCAM-1 samples show linear curves that were parallel to the standard curves obtained using the kit standards. These results indicate that the immunoassay kit can be used to determine relative mass values for natural human NCAM-1.

**PRINCIPLE OF THE ASSAY**

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for NCAM-1 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any NCAM-1 present is bound by the immobilized antibody. After washing away any unbound substances, a biotinylated polyclonal antibody specific for NCAM-1 is added to the wells. Following a wash to remove any unbound antibody-biotin reagent, HRP link Streptavidin is added to the wells. After washing away any unbound enzyme, a substrate solution is added to the wells and color develops in proportion to the amount of NCAM-1 bound in the initial step. The color development is stopped and the intensity of the color is measured.

**LIMITATIONS OF THE PROCEDURE**

- \_ FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- \_ The kit should not be used beyond the expiration date on the kit label.
- \_ Do not mix or substitute reagents with those from other lots or sources.
- \_ It is important that the dilution buffer selected for the standard curve be consistent with the samples being assayed.
- \_ If samples generate values higher than the highest standard, dilute the samples with Dilution Buffer and repeat the assay.
- \_ Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- \_ This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors

have been tested in the immunoassay, the possibility of interference cannot be excluded.

**MATERIALS PROVIDED**

Description	Code	Quantity
<b>Microplate</b> - 96 well polystyrene microplate (12 strips of 8 wells) coated with a mouse monoclonal antibody against NCAM-1.	266-01-01	1 plate
<b>NCAM-1 Standard</b> – 4000 pg/vial of recombinant human NCAM-1 in a buffered protein base with preservatives; lyophilized.	266-01-02	1 vial
<b>Detection Antibody Concentrate</b> – 140 µL / vial, 100-fold concentrated of Biotinylated polyclonal antibody against NCAM-1 with preservatives; lyophilized.	266-01-03	1 vial
<b>Positive Control</b> - one vial of recombinant human NCAM-1, lyophilized	266-01-04	1 vial
<b>Streptavidin-HRP Conjugate</b> - 75 ul/vial, 240-fold concentrated solution of Streptavidin conjugate to HRP	SAHRP	1 vial
<b>Dilution Buffer</b> - 60mL/vial of buffered protein based solution with preservatives	DB01	1 vial
<b>Wash Buffer</b> -50 ml/vial, 10-fold concentrated buffered surfactant, with preservative.	WB01	1 vial
<b>TMB Substrate Solution</b> - 11ml/vial of TMB substrate solution	TMB01	1 vial
<b>Stop Solution</b> (0.5M HCl) , 13ml/vial of 0.5M HCl	S-STOP	1 vial
<b>Plate Covers</b> – Plate sealer	EAPS	1

**STORAGE**

**Unopened Kit:** Store at 2 - 8° C for up to 5 months. For longer storage, unopened Standard, Positive Control and Detection Antibody Concentrated as

well as Dilution Buffer BD09 should be stored at -20 or -70 °C. Do not use past kit expiration date.

**Opened / Reconstituted Reagents:** Reconstituted Standard , Antibody Solution SHOULD BE STORED at -20 °C or – 70°C for up to one months. Anti Rabbit IgG - HRP Conjugate 100-fold concentrated and other components may be stored at 2 - 8°C for up to 6 months.

**Microplate Wells:** Return unused wells to the plastic bag containing the desiccant pack, reseal along entire edge of zip-seal. May be stored for up to 6 months at 2 - 8° C.

**OTHER SUPPLIES REQUIRED**

- Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 540 nm or 570 nm.
- Microplate shaker (250-300rpm).
- Pipettes and pipette tips.
- Deionized or distilled water.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- 100 mL and 500 mL graduated cylinders.

**PRECAUTIONS FOR USE**

All reagents should be considered as potentially hazardous. The stop solution contains diluted Hydrochloric acid. Appropriate care, therefore, should be taken while handling this solution. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water.

**SAMPLE COLLECTION AND STORAGE**

**Plasma** - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤-20° C. Avoid repeated freeze-thaw cycles. Optimal dilutions should be determined by each laboratory for each application.

**REAGENT PREPARATION**


**Bring all reagents to room temperature before use.**

**Wash Buffer** - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved.

Dilute 50 mL of Wash Buffer Concentrate into deionized or distilled water (450 mL) to prepare 500 mL of Wash Buffer.

**NCAM-1 Standard - Refer to vial label for reconstitution volume.** Reconstitute the **NCAM-1 Standard** with 1 ml of Dilution Buffer. This reconstitution produces a stock solution of 4000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 125µL of Dilution Buffer into tube #1, then transfer 375µL of 4000 pg/mL stock solution to tube #1. This will create the high standard of 3000pg/mL). Then, pipette 250µL of Dilution Buffer into tubes #2 to #7. Use the high standard solution to produce a dilution series (see below). Mix each tube thoroughly before the next transfer. The 3000 pg/mL standard serves as the high standard. The Dilution Buffer serves as the zero standard (0 pg/mL).

Standard	Standard	Reagent Diluent	Concentration
Stock	Powder	1000 µl	4000 pg/mL
# 1	375 µl of stock	125 µl	3000 pg/mL
# 2	250 µl of 1	250 µl	1500 pg/mL
# 3	250 µl of 2	250 µl	750 pg/mL
# 4	250 µl of 3	250 µl	375 pg/mL
# 5	250 µl of 4	250 µl	187.5 pg/mL
# 6	250 µl of 5	250 µl	93.75 pg/mL
#7	250 µl of 6	250 µl	46.875 pg/mL

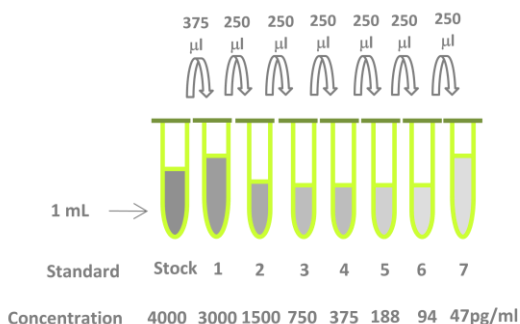

**Positive Control** - Reconstitute the **Positive Control** with 1.0 mL Dilution Buffer. In a micro-centrifuge tube, add 100 µL Dilution Buffer then transfer 300µL of the stock Positive Control solution to the tube to prepare the Positive Control working solution. *Positive Control should be prepared and used immediately.*

**Detection Antibody-** Reconstitute the **Detection Antibody Concentrate** with 140 µL of Dilution Buffer to produce a 100-fold concentrated stock solution. Pipette 13.86 mL of Dilution Buffer into a 15 ml centrifuge tube and transfer 140 µL of 100-fold concentrated stock solution to prepare working solution.

**Streptavidin-HRP Conjugate** - Pipette 12 mL of Dilution Buffer into a 15 ml centrifuge tube and transfer 50 µL of the 240-fold concentrated stock solution to prepare working solution.

**ASSAY PROCEDURE**

**Bring all reagents and samples to room temperature before use. It is recommended that standards be assayed in duplicates.**


1. Prepare all reagents and working standards as directed in the previous sections.
2. Remove excess micro-plate strips from the plate frame, return them to the foil pouch containing the desiccant pack, reseal.
3. Add 100 µL of Dilution Buffer to Blank well (A2,A3).
4. Add 100 µL of Standard (from B2,B3 to G2,G3 and G4,G5), sample, or control (F4,F5) per well. Cover with plate sealer. Incubate for 2 hours on micro-plate shaker at room temperature. A plate layout is provided to record standards and samples assayed.
5. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (300 µL) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
6. Add 100 µL of Detection Antibody working solution to each well. Cover with plate sealer. Incubate for 2 hours on micro-plate shaker at room temperature.
7. Repeat the aspiration/wash as in step 5.
8. Add 100 µL of **Streptavidin-HRP Conjugate** working solution to each well. Incubate for 45 minutes on micro-plate shaker at room temperature.
9. Repeat the aspiration/wash as in step 5.
10. Add 100 µL of Substrate Solution to each well. Incubate for 5-15 minutes at room temperature. **Protect from light.**

11. Add 100 µL of Stop Solution to each well. The color in the wells should change from blue to yellow. If the color in the wells is green or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
12. Determine the optical density of each well within 15 minutes, using a micro-plate reader set to 450 nm.

**CALCULATION OF RESULTS**

Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density. Create a standard curve by reducing the data using computer software capable of generating a log-log curve fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the NCAM-1 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

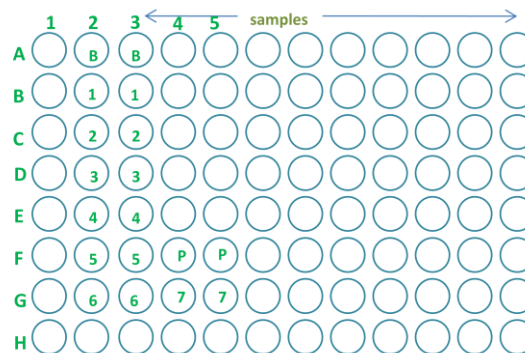

**TYPICAL DATA**

These standard curves\* are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

NCAM-1(pg/ml)	Average OD450 (Corrected)*
46.875	0.028
93.75	0.065
187.5	0.131
375	0.259

750	0.493
1500	0.993
3000	1.983

\*Lot No.:

\*\* Positive Control: 331-554 pg/ml

**CALIBRATION**

This immunoassay is calibrated against a highly purified NS0-expressed recombinant human NCAM-1.

**SENSITIVITY**

Twenty-five assays were evaluated and the minimum detectable dose (MDD) of NCAM-1 was 24.5 pg/mL.

**SPECIFICITY**

This assay recognizes both natural and recombinant human NCAM-1. The factors listed below were prepared at 50 ng/mL in Dilution Buffer, and assayed for cross reactivity. Preparations of the following factors at 50 ng/mL in a mid-range rh NCAM-1 control were assayed for interference. No significant cross-reactivity or interference was observed.

Human Recombinant Proteins:

NCAM-L1/Fc Chimera

BCAM/Fc Chimera

EpCAM/TROP-1/Fc Chimera

ALCAM-L1/Fc Chimera

Mouse Recombinant Proteins:

MadCAM-1/Fc Chimera

**SUMMARY OF ASSAY PROCEDURE**

Prepare reagents, samples and standards
↓
Add 100µl of standard, samples, positive control to each well. Incubate 2 hours on the plate shaker at RT.
↓
Aspirate and wash 4 times.
↓
Add 100 µl Detection Antibody to each well. Incubate 2 hours on the plate shaker at RT.
↓
Aspirate and wash 4 times.
↓
Add 100 µl Streptatvin HRP conjugate to each well. Incubate 45 minutes on the plate shaker at RT.
↓
Aspirate and wash 4 times.
↓
Add 100 µl Substrate to each well. Incubate 5-15min on the bench top. Protect from light.
↓
Add 100 µl Stop Solution to each well. Read 450nm within 15 min